




Open Linear Stapler



Available in four different closing lengths - 30mm, 45mm, 60mm and 90mm.

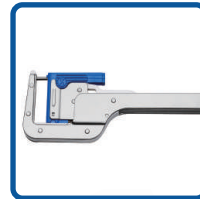


The locating pin can be placed manually or automatically, providing doctors with a variety of options for easy surgical operation.




Lockable handle - the handle locks in the back position after firing which is helpful for the surgeon to determine if the device is fired, hence avoiding miss-firing.



Two different staple heights: 3.5mm and 4.8mm available, suitable for different tissues.

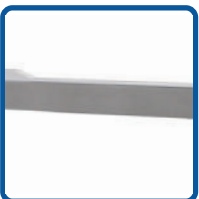

Long narrow shaft which allows for greater access and visibility.

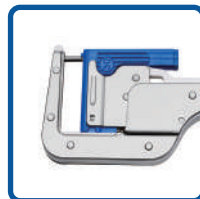

When cutting tissue, the cutting guide can help ensure safe and convenient cut.

Unique safety feature!

The colour indicator window shows different colours for different firing phases and reduces the possibility of misoperation.



Primary status
 semi-closed status
 fully-closed status
 post-firing status



Product Codes

Open Linear Stapler



Product Code	Description	Staple Length	Staple Height
FLSLF30-3.5	Linear Stapler	30mm	3.5mm
FLSLF30-4.8	Linear Stapler	30mm	4.8mm
FLSLF45-3.5	Linear Stapler	45mm	3.5mm
FLSLF45-4.8	Linear Stapler	45mm	4.8mm
FLSLF60-3.5	Linear Stapler	60mm	3.5mm
FLSLF60-4.8	Linear Stapler	60mm	4.8mm
FLSLF90-3.5	Linear Stapler	90mm	3.5mm
FLSLF90-4.8	Linear Stapler	90mm	4.8mm

Cartridges for Linear Cutter Stapler



Product Code	Staple Length	Staple Height
FLSLH3.5A	30mm	3.5mm
FLSLH4.8A	30mm	4.8mm
FLSLH3.5B	45mm	3.5mm
FLSLH4.8B	45mm	4.8mm
FLSLH3.5C	60mm	3.5mm
FLSLH4.8C	60mm	4.8mm
FLSLH3.5E	90mm	3.5mm
FLSLH4.8E	90mm	4.8mm